

USA E8 CENTRIFUGES



E8 TOUCH



E8 VARIABLE



E8 FIXED

ADVANTAGES

Accurate speeds for specific g-forces
Proper separations without cell damage
Quiet, vibration-free operation
Auto-brake stops the rotor in less than 30 seconds
Designed, Engineered, and Assembled
Proudly in America

FEATURES

Angled rotor, 8-tube capacity
3,500 rpm
Holds 3ml to 15ml size tubes
12 volt DC maintenance-free motor
Clear view port in lid for using tachometer
Suction-cupped feet prevent slipping
110/220v auto-switching power adapter
(CE, UL, cUL approved adapter)
Lid safety shut-off switch
Maximum RCF=1,534g
E8 Fixed & E8 Variable: Auto-off 30-minute timer with bell
E8 Touch: Touch screen, programmable

WEIGHT AND DIMS

Height: 9.5" (241mm)
Length: 11" (279mm)
Width: 11" (279mm)
Weight: 9.5 lbs (4.3kg)



The USA E8 Centrifuge is designed, engineered, and assembled at LW Scientific's ISO-certified, FDA-registered manufacturing facility in Lawrenceville, GA. Expect American quality as well as customary LW Scientific support. This model is loaded with features, such as higher speeds for faster and cleaner separations, quieter operation for peaceful work environments, lower chamber temperatures for cooler test-tube samples, shorter height to fit under any cabinet, and auto-calibrating speed function for precision results. It even has an auto-brake to stop the rotor in less than 30 seconds to save time in busy labs.



Designed, Engineered &
Assembled Proudly in America

STANDARD E8

INCLUDES

Eight-place rotor
Eight 15ml tube sleeves
Eight 13x75mm tube sleeve inserts




F	E8C-U8AF-1503	USA E8 Fixed-Speed, Angled 8-Place Centrifuge with Timer
V	E8C-U8AV-1503	USA E8 Variable-Speed, Angled 8-Place Centrifuge with Timer
T	E8C-U8AT-15T3	USA E8 Digital-Speed, Touch Screen, Angled 8-Place Centrifuge

COMBINATION E8

INCLUDES

Eight-place rotor
Eight 15ml tube sleeves
Eight 13x75mm tube sleeve inserts
Six green tube inserts to spin 1ml-2ml
Two 4-place Crit Carrier inserts
(Spins 75mm capillary tubes)
EZ Reader Card

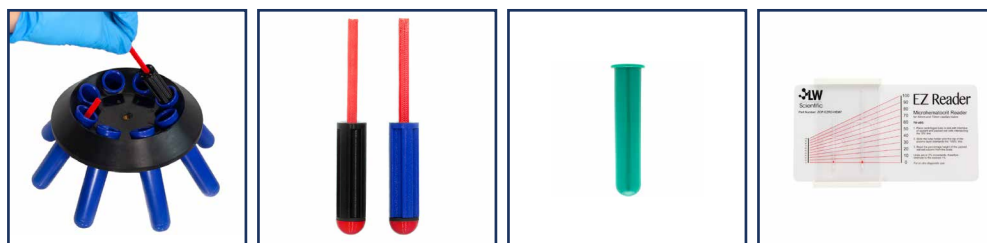

F	E8C-U8AF-CR03	USA E8 Fixed-Speed, Angled 8-Place Centrifuge with Crit Carriers
V	E8C-U8AV-CR03	USA E8 Variable-Speed, Angled 8-Place Centrifuge with Crit Carriers
T	E8C-U8AT-CRT3	USA E8 Digital-Speed, Touch Screen, Angled 8-Place Centrifuge with Crit Carriers

PORTABLE E8

INCLUDES

Eight-place rotor
Eight 15ml tube sleeves
Eight 13x75mm tube sleeve inserts
Car seat
Car adapter cord



F	E8C-U8AF-150P	USA E8 Fixed-Speed, Angled 8-Place Centrifuge, Portable
V	E8C-U8AV-150P	USA E8 Variable-Speed, Angled 8-Place Centrifuge, Portable
T	E8C-U8AT-15TP	USA E8 Digital-Speed, Touch Screen, Angled 8-Place Centrifuge, Portable